

dragonfly[®] crystal: an ideal solution to problematic liquids in screening environments

Joby Jenkins, Gillian Lewis
SPT Labtech, Melbourn Science Park, Melbourn, Royston, Herts, SG8 6HB

introduction

The accurate pipetting of viscous liquids and cross-contamination between wells can be highly problematic when using conventional pipetting technologies, such as air displacement pipetting.

The dragonfly[®] crystal positive displacement liquid handling instrument eliminates these issues and can be used to optimise assays and formulations in a wide variety of application areas including drug discovery, biotechnology, clinical and food industries.

Very high or very low viscous solutions, such as 100% glycerol, PEGs, alcohols, or detergent solutions, including thixotropic non-Newtonian ones, are easily pipetted with very high accuracy and precision. Liquids are dispensed via a non-contact unique dispensing technology to ensure there is zero cross-contamination between wells and the disposable tips eliminate cross-contamination between samples.

1. dragonfly crystal: digital liquid handler

The novel dispensing technology of the dragonflycrystal digital liquid handler enables dispensing directly into assay plates with a very large dynamic range of 0.5 μ L to 4 mL. Independent volume control and simultaneous digital dispensing from up to 10 pipetting tips provide rapid dispensing; less than 5 minutes for dispensing 4 solutions into a 96-well plate and total freedom of plate layout. Suitable applications range from formatting multi-dimensional reagents, optimisation experiments to backfilling of each well with assay buffer or PCR master mix.



2. pipetting technology

A dragonfly crystal is equipped with an array of either 5 or 10 pipetting heads. These utilise positive displacement pipettes providing reliable and accurate dispensing of a broad range of liquids without any classification (Fig 1: A).

The use of disposable syringe and plunger guarantees no cross-contamination. Each pipette works independently to dispense any volume into any well for complete assay flexibility. Each liquid is aspirated from a dedicated easy-to-fill reservoir (Fig 1: B & C).

Unused liquid at the end of a run can thus be recovered resulting in a dead volume of < 0.5 μ L. The whole process is achieved without valves, air gaps or system fluid resulting in high reliability.


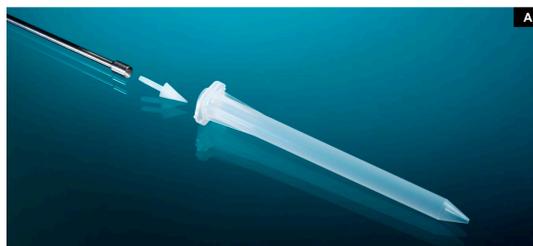


Fig 1: dragonfly crystal's disposable positive displacement pipettes for non-contact dispensing

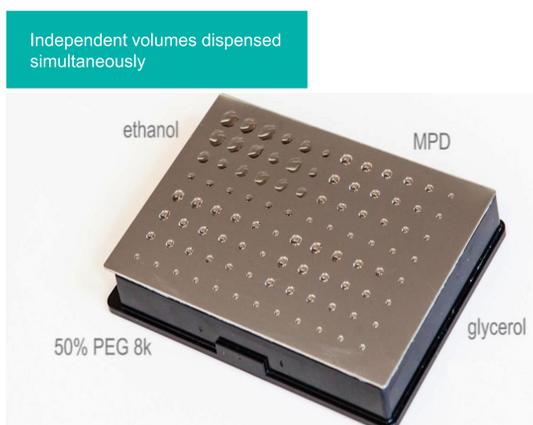

3. result of pipetting 100% glycerol: a Newtonian liquid

Fluids may be classified as Newtonian or non-Newtonian. In a Newtonian fluid, there is a linear relation between the magnitude of applied shear stress and the resulting rate of deformation. Newtonian fluids have a constant viscosity (dynamic or absolute viscosity) at a given temperature such as water, ethanol or aqueous solutions of salts and sugar.

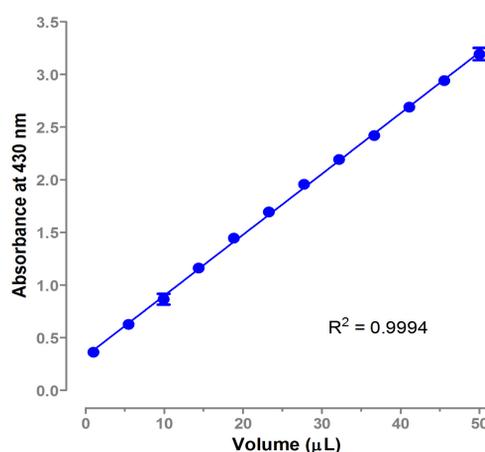


Fig 2: Linearity of dispense for 100% glycerol

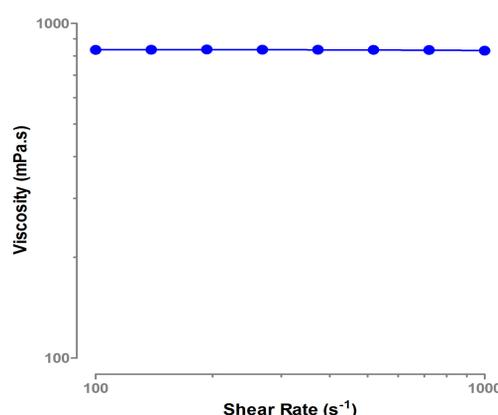


Fig 3: Viscosity vs shear rate of 100% glycerol. The method we used was a rotational rheometer with a cone and plate.

4. result of pipetting a non-Newtonian liquid

A wide range of industrially important liquids, such as solutions of high molecular weight polymers, colloids, suspensions, and emulsions exhibit more complex behaviour, which is termed non-Newtonian. In non-Newtonian fluid, there is a nonlinear relationship between the magnitude of applied shear stress and the rate of angular deformation.

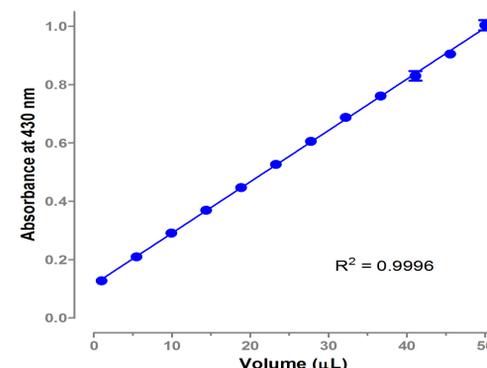


Fig 4: Linearity of dispense of 17.5 mM dyed household detergent using dragonfly crystal

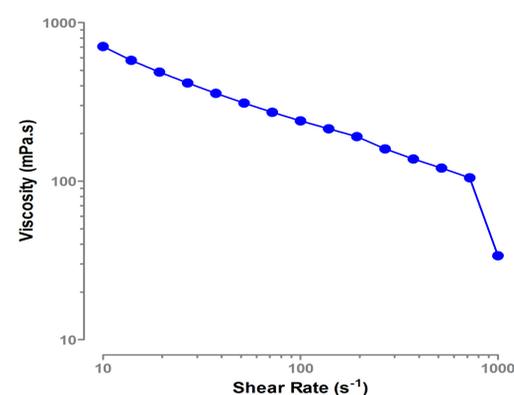


Fig 5: Viscosity vs shear rate. The method we used was a rotational rheometer with a cone and plate

5. pipetting data-accuracy and precision

liquid type	% inaccuracy	% CV
4 μ L glycerol	0.18	1.20
4 μ L Lutensit	2.50	0.80
4 μ L Pluriol P900	2.0	0.73
1 μ L Tween 80	1.30	0.91
1 μ L Tween 20	1.20	0.89



conclusions

Our work demonstrates:

- dragonfly crystal's ability to pipette difficult liquids that often cause blocking and contamination in an optimisation study
- the dragonfly crystal liquid handler is valuable to the discovery workbench, eliminating the tedium of complicated microplate set-ups while maintaining flexibility